




BIONIT LABS®

TURNING DISABILITIES INTO NEW POSSIBILITIES



ADAM'S HAND®

The First of its Kind

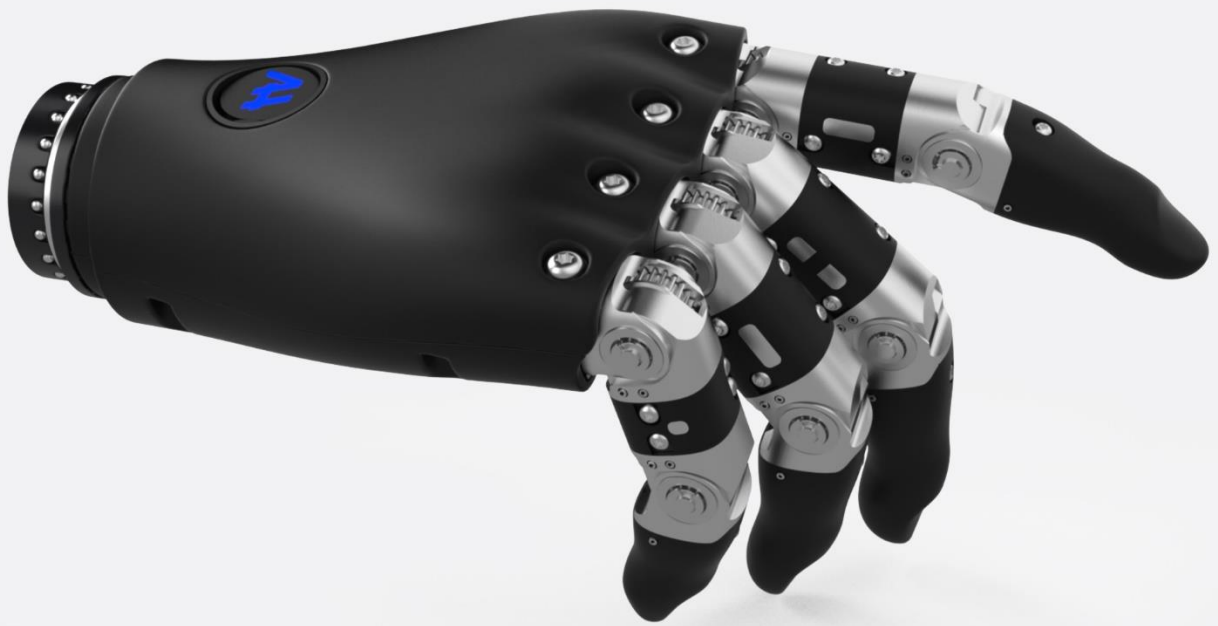


ADAM'S HAND
Mod. AH02-CDY-EE

DATASHEET



WARNING: the image shown on the cover contains the prosthetic hand without a glove, but for representative purposes only, to highlight the actual position of the button and the quick disconnect wrist. For the safety of the user and for the correct functionality of the "Adam's Hand®" prosthetic hand, this must ALWAYS be used with its own "Adam's Skin" glove mod. AS02-CDY-EE.

COMMERCIAL NAME	Adam's Hand®	
REF	AH02-CDY-EE, with C = side code, D = size code, Y = wrist code, EE = colour code	
INTENDED USE	The "Adam's Hand®" prosthetic hand is a non-invasive medical device intended to be used in myoelectric prostheses made for subjects with mono- or bi-lateral amputations, starting from a transradial amputation level.	
TECHNICAL FEATURES	COMPATIBILITY	
	Consult the chart on p. 4, containing the list of devices compatible with the Adam's Hand®.	
	PRODUCT DATA	
	Size	Medium (M)
	Side	Left (C=L) Right (C=R)
	N° actuators	2
	Wrist type	Quick Disconnect Wrist (Y=1)
	Weight	560 g
	Thumb adduction/abduction	Passive, 3 selectable positions: 0° (lateral), 45° (neutral), 90° (opposition).
	Expected lifetime	6 years
	Warranty	3 years (with possibility of extension up to 5 years)
	Available Colours	Black (EE=B1), White/Transparent (EE=W1) Skin tone 1 (EE=N1) Skin tone 2 (EE=N2) Skin tone 3 (EE=N3) Skin tone 4 (EE=N4) Skin tone 5 (EE=N5) Skin tone 6 (EE=N6) Skin tone 7 (EE=N7) Skin tone 8 (EE=N8)
	ELECTRICAL AND MECHANICAL FEATURES	
	Supply voltage	6 ÷ 8.4 V (DC)
	Supply current	6 A
Degree of protection against liquid and solid particles penetration (IEC 60529)	IP67	
Max. grip strength	<ul style="list-style-type: none"> Lateral grip: 75 N Tridigital grip: 50 N Power grip: 100 N 	

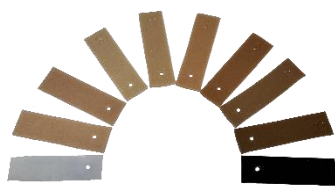

FUNCTIONALITIES	Closing time	<ul style="list-style-type: none"> Lateral grip: 0.75 s Tridigital grip: 0.7 s Power grip: 1.2 s
	Vertical push-down load through knuckles	100 kg
	Max. static load	<ul style="list-style-type: none"> With closed hand (hook grip): 35 kg On a single finger: 16 kg
	Static force on the chassis	75 kg
	SOFTWARE	
	<p>The Adam's Hand® prosthetic hand integrates a Bluetooth Low Energy (BLE)® module, useful for connection to an external device hosting the “MyoLogic” (iOS and Android app) or “MyoLogic Pro” (iOS app) Apps, respectively dedicated to Patient and Orthopedic Technician or authorized Distributor, through which it is possible to read the electromyographic signals detected by the electrodes, calibrate the operating parameters and calibrate the Adam's Hand® to adapt it to the individual's muscle tone. Once the Adam's Hand® has been installed and calibrated, the use of these apps for the use of the device by the Patient and the Orthopedic Technician or Authorized Distributor is recommended but not required.</p>	
	Minimum system requirements	MyoLogic: iOS® 15 or next MyoLogic Pro: iPadOS® 13 or next
	OPERATING CONDITIONS	
	Use temperature	from -10 °C/+14 °F to + 40 °C/+104 °F
	Usage humidity	30% ÷ 85 %, non-condensing
Storage and transport temperature	from -20 °C/-4 °F to +70 °C/+158 °F	
Relative humidity of storage and transport	≤ 85 %, non-condensing	
MULTIFUNCTION BUTTON		
<p>By activating the button on the back of the hand, the User can perform the following functions, which can be programmed via the App:</p> <ul style="list-style-type: none"> activation/deactivation of the EMERGENCY mode activation/deactivation of the STANDBY mode modification of the Mode of Use (ECO, NORMAL, BOOST) activation/deactivation of the Bluetooth Low Energy connection activation/deactivation of the Proportional Control new Calibration of the device. 		
VISUAL ALERTS		
<p>Polychromatic LED integrated in the button on the back of the hand for signaling alarms, device status changes and confirmation of activations/deactivations implemented with the button or via the App.</p>		
ACOUSTIC ALERTS		
<p>4 kHz frequency buzzer for signalling alarms, device status changes and confirmation of activations/deactivations implemented with the multifunction button or via the App.</p>		
HAPTIC ALERTS		
<p>The vibration is proportional to the grip force exerted by the prosthetic hand. It is also used for signalling alarms, device status changes and confirmation of the activation/deactivation of the prosthetic hand functions (implemented with the multifunction button or via the App).</p>		


	SENSORS
	<ul style="list-style-type: none"> • Internal temperature sensors • Automatic protection against: overtemperature, overcurrent and overvoltage • Sending data via Bluetooth Low Energy • The inertial platform (IMU) integrated in the Adam's Hand® allows the device to have information on the spatial attitude (without measurement function) through acceleration and rotation sensors with 16-bit resolution.
TESTS	Each medical device is tested before the shipment, in accordance with company procedures. The reference standards are affixed on the declaration of conformity attached to the device.
TECHNICAL RULES	IEC 60601-1, IEC 60601-1-2, IEC 60601-1-11, IEC 62304, IEC 62366 ISO 10993-10, ISO 22523 RED ETSI 300 328
LABELLING	Labelling in accordance with UNI EN ISO 15223-1, IEC 60601-1; copy of the label is present in the Technical Manual and User Manual.
DISPOSAL INSTRUCTIONS	 <p>This medical device must be managed in accordance with art. 13 - Legislative Decree 25 July 2005, n. 151 "Implementation of directives 2002/95/CE, 2002/96/CE and 2003/108/CE, relating to the reduction of the use of dangerous substances in electrical and electronic equipment, as well as waste disposal".</p>

Chart 1: Devices compatible with the Adam's Hand ® prosthetic hand.

Manufacturer	Product	Type	Model / Part Number
BionIT Labs® S.r.l.	ThunderCell Battery	External battery	EB01 EB02-E
BionIT Labs® S.r.l.	Ricettacolo con Cavo per ThunderCell Battery	Receptacle for external battery	ER01-E-XX
BionIT Labs® S.r.l.	Wave Electrode	Analogic electrode	AE02-NN
Otto Bock HealthCare GmbH	Coupling piece	Wrist interface	10S4
Otto Bock HealthCare GmbH	Coaxial plug (4 bands)	Wrist interface	9E169
Otto Bock HealthCare GmbH	Lamination Ring	Wrist interface	10S1
Steeper Group	Quick Disconnect Wrist Unit	Wrist interface	QDAWA-40 / -45 / -50
Motion Control	Coaxial plug (4 band)	Wrist interface	1701009
Motion Control	Endcap	Wrist interface	3010758
Motion Control	Lamination Collar for MC Wrist	Wrist interface	1100288
Össur	Coupling piece	Wrist interface	PL091032
Össur	4-way coaxial plug	Wrist interface	PL091036
Össur	Lamination Ring	Wrist interface	091037
Otto Bock HealthCare GmbH	Energy Pack®	External battery	757B20
Otto Bock HealthCare GmbH	Energy Pack®	External battery	757B21
Otto Bock HealthCare GmbH	MyoEnergy Integral®	Internal battery	757B35=0 / =1 / =3 / =4 / =5
Otto Bock HealthCare GmbH	MyoBock® Electrode	Analogic electrode	13E200=50 / =60
Otto Bock HealthCare GmbH	Suction Socket Electrode	Analogic electrode	13E202=50 / =60
Össur	1300/2000 mAh Battery	Internal battery	PL000335/6
Össur	Rechargeable Li-Polymer Battery	Internal battery	000172B
Össur	Replaceable Battery	External battery	SA069313
Össur	Compact Electrode Kit	Analogic electrode	PL091-XXX
Steeper Group	Lithium Polymer Battery	External battery	BLPA72
Steeper Group	S-Charge System®	Internal battery	SC2200
Steeper Group	Electrode	Analogic electrode	ELEC50 / ELEC60
TASKA	TASKA® Battery system	Internal battery	-
COVVI Ltd.	Power Supply systems	Internal battery	COVSB-1600
COVVI Ltd.	Power Supply systems	Internal battery	COVSB-2600
COVVI Ltd.	Electrode	Analogic electrode	CEL-50 / CEL-60
IBT	FlexCell®	Internal battery	2017200
IBT	FlexCell Mini®	Internal battery	2037200
Vincent Systems GmbH	VINCENTaccu flex®	Internal battery	-

Manufacturer's contacts

For any information, request or complaint, please contact:



BionIT Labs S.r.l.

Via Cracovia, 1

73010 – Soletto (LE) – ITALY

Certified Company





BIONIT LABS[®]
TURNING DISABILITIES INTO NEW POSSIBILITIES



+39 0836 307525



www.bionitlabs.com



support@bionitlabs.com



Via Cracovia, 1
73010 Soleto (LE) - Italia